# ToothWave

H7001

- EN USER MANUAL
- DE BEDIENUNGSANLEITUNG
- FR MODE D'EMPLOI
- NL GEBRUIKSAANWIJZING
- ES MANUAL DEL USUARIO
- PT MANUAL DO UTILIZADOR
- GR ΕΓΧΕΙΡΙΔΙΟ ΧΡΗΣΤΗ
- TR KULLANIM KILAVUZU
- IT MANUALE PER L'UTENTE

Patent protected  
technology by

TOOTH  
WAVE™

SCIENCE DRIVEN ORAL CARE

DentalRF™  
Technology

Silk'n®


# Silk'n ToothWave


#10

#9
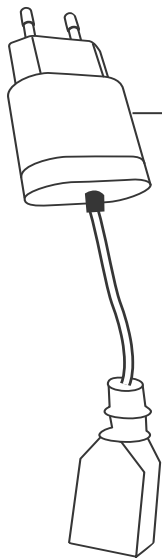
#8
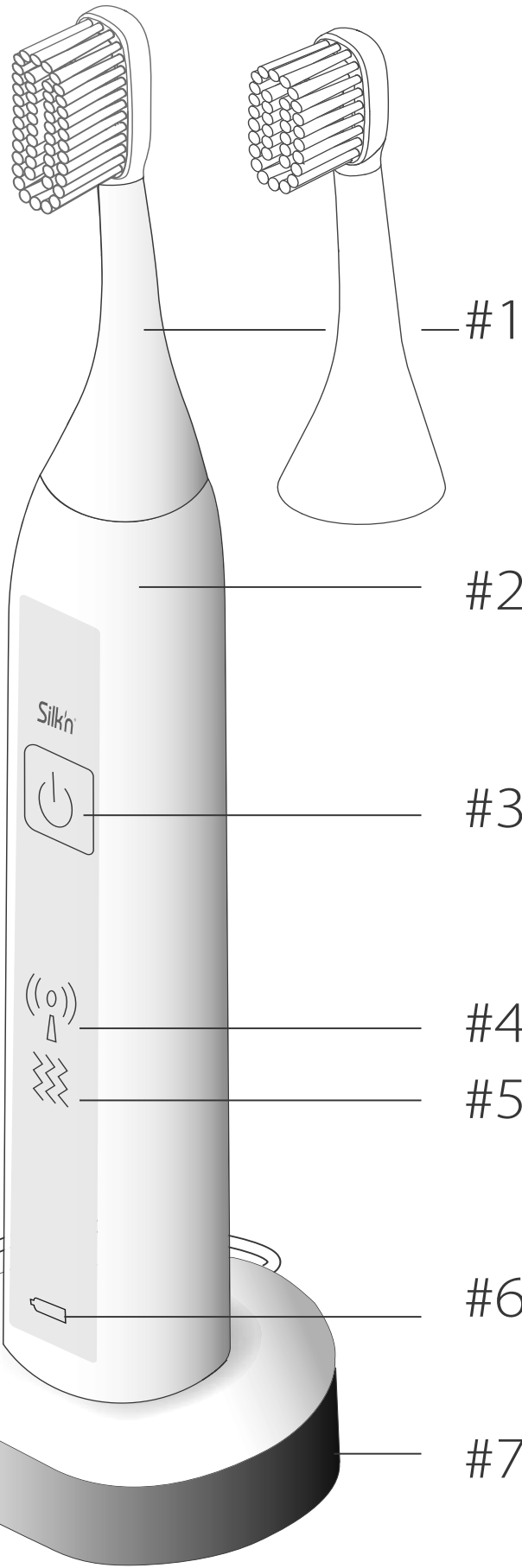
# 1. Warnings & contraindications

---

---

EN

A statement or event that indicates a possible or impending danger, problem, or other unpleasant situation when using the device.

## Electricity & safety warnings

- Always ensure that your hands are dry before connecting the charging cradle to the electricity supply.
- Do not place or store the charging cradle and the USB wall adapter where it can fall, or be pulled into any liquid. The charging cradle and the USB wall adapter are not waterproof.
- For the purpose of recharging the battery, only use the USB wall adapter provided with this device.
- Do not use the charging cradle or the USB wall adapter if either has fallen into liquid.
- If the charging cradle malfunctions, unplug immediately.
- Never use the charging cradle or the USB wall adapter if they are damaged or aren't working properly.

- Keep the device, all parts and accessories away from heated surfaces.
- Do not attempt to modify, open or repair the device. You may expose yourself to dangerous electrical components and may cause serious harm to yourself. It will also void your warranty.
- Do not use if the parts or accessories are not supplied by the manufacturer as defined in this manual, or if they are damaged in any way, do not appear to be working properly, or if you see or smell smoke. In such cases, stop using the device, and contact customer service.
- The manufacturer accepts no liability for damage or injuries caused by improper or incorrect use.
- Do not clean any of the parts in the dishwasher.
- Children shall not play with the device. Cleaning and user maintenance shall not be made by children.
- Children must not play with the packaging. They may injure themselves with the packaging material, or swallow small parts and choke.

- This device contains batteries that are non-replaceable.
- Do not dispose of the device or the battery with household waste! The law requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices.
- Dispose of the battery safely.
- If the product reaches its end of life, the battery must be removed before disposal (see section 12. Battery removal).

## **Contraindications**

This device including all its parts must not be used by the following people: children below the age of 18, adults with limited physical, sensory or psychological capacities, those lacking experience or knowledge in how to use the device in a safe way, or those who do not understand the hazards involved.

### **Do not use if you:**

- have a pacemaker, internal defibrillator, another active implanted device, or if you have any medical concerns about using this device.

- have a current or past history of oral cancer or oropharyngeal cancer, or if you have any other type of cancer, or pre-malignant lesions.
- have an impaired immune system due to immuno-suppressive diseases such as HIV, or you use immuno-suppressive medication.

### **Consult your physician before use if you:**

- have concurrent conditions, such as cardiac disorders, epilepsy, uncontrolled hypertension, and liver or kidney diseases.
- have poorly controlled endocrine disorders, such as diabetes or thyroid abnormality.
- are pregnant or nursing.

### **Consult your dentist before use if you:**

- have had oral or gum surgery in the last two months.
- have severely injured gums.

### **Additional warnings:**

- Do not bite the brush head.

## 2. Signs & symbols

---

EN

	Important information about using or maintaining the device.
	Products marked with this symbol fulfill the requirements of the EU directives.
IP67	Protected from immersion in water with a depth of up to 1 meter (or 3.3 feet) for up to 30 mins.
	Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE Directive.
 Li-Ion	Lithium inside.
	The built-in rechargeable ToothWave battery contains substances that may pollute the environment. Always remove the battery before you discard at an official collection point for batteries.

## 3. Get to know Silk'n ToothWave

---

### 1.1 Intended use

Silk'n ToothWave is intended for the removal of teeth discolouration, stains, plaque and tartar (calculus). Over time and repeated use, the reduction of tartar causes the remission of gum disease (gingivitis) and bleeding, and the prevention of tooth decay. See our websites and videos for more information about our patented DentalRF™ technology.

Silk'n ToothWave is for personal home-use only, and is not intended for use on multiple patients in a dental practice or institution. Do not use for any purpose other than as described in this manual.

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© 2019 Home Skinovations Ltd. Minden jog fenntartva.